College of Engineering & Technology, SRMIST, Kattankulathur – 603203



Program: B. Tech. (CSE, IT, SWE and Nano Technology)

18PYB103J - Semiconductor Physics - Instructional Laboratory Manual

EX. NO. NAME OF THE EXPERIMENTS PAGE NO.

1.	DETERMINATION OF HALL COEFFICIENT AND CARRIER TYPE FOR A SEMI-CONDUCTING MATERIAL	3
2.	BAND GAP DETERMINATION USING POST OFFICE BOX	6
3.	TO STUDY V-I CHARACTERISTICS OF A LIGHT DEPENDENT RESISTOR (LDR)	9
4.	RESISTIVITY DETERMINATION FOR A SEMICONDUCTOR WAFER USING FOUR PROBE METHOD	12
5.	STUDY OF V-I AND V-R CHARACTERISTICS OF A SOLAR CELL	15
6.	CHARACTERISTIC OF PN JUNCTION DIODE UNDER FORWARDBIAS	17
7.	TO STUDY ILLUMINATION AND V-I CHARACTERISTICS OF A PHOTO CELL	19
8.	PARTICLE SIZE DETERMINATION USING LASER	22
9.	STUDY OF ATTENUATION AND PROPAGATION CHARACTERISTICS OF OPTICAL FIBER CABLE	24
10.	CHARACTERISTIC OF PN JUNCTION DIODE UNDER REVERSE BIAS	28
11.	DETERMINATION OF EFFICIENCY OF SOLAR CELL	30
12.	CALCULATION OF LATTICE CELL PARAMETERS – X-RAY DIFFRACTION	33
13.	MINI PROJECT – CONCEPT BASED DEMONSTRATION	36

1. Determination of Hall Coefficient and carrier type for a Semi-conducting Material

Aim

To determine the hall coefficient of the given n type or p-type semiconductor

Apparatus Required

Hall probe (n type or p type), Hall effect setup, Electromagnet, constant current power supply, gauss meter etc.,

Principle

Hall effect: When a current carrying conductor is placed in a transverse magnetic field, a potential difference is developed across the conductor in a direction perpendicular to both the current and the magnetic field.

Formula

i) Hall coefficient (R_H) = $\frac{V_H.t}{IH} \times 10^8 \text{ cm}^3 \text{ C}^{-1}$

where V_H = Hall voltage (volt)

t = Thickness of the sample (cm)

I = Current (ampere)

H = Magnetic field (Gauss)

ii) Carrier density (n)= $\frac{1}{R_H q}$ carriers / cm⁻³

where $R_H = Hall coefficient (cm^3 C^{-1})$

q = Charge of the electron or hole (C)

iii) Carrier mobility (μ)= $R_H \sigma \text{ cm}^2 \text{V}^{-1} \text{ s}^{-1}$

where $R_H = \text{Hall coefficient } (\text{cm}^3 \text{C}^{-1})$

 σ = Conductivity (C V⁻¹s⁻¹cm⁻¹)

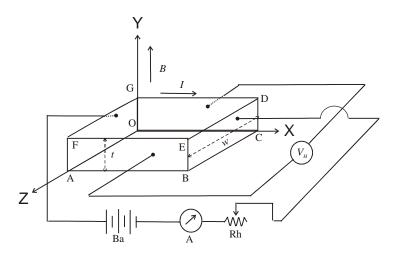


Fig.1.1 Hall Effect Setup

Procedure

- 1. Connect the widthwise contacts of the hall probe to the terminals marked as 'voltage' (i.e. potential difference should be measured along the width) and lengthwise contacts to the terminals marked (i.e. current should be measured along the length) as shown in fig 1.1.
- 2. Switch on the Hall Effect setup and adjust the current say 0.2 mA.
- 3. Switch over the display in the Hall Effect setup to the voltage side.
- 4. Now place the probe in the magnetic field as shown in fig and switch on the electromagnetic power supply and adjust the current to any desired value.Rotate the Hall probe until it become perpendicular to magnetic field. Hall voltage will be maximum in this adjustment.
- 5. Measure the hall voltage and tabulate the readings.
- 6. Measure the Hall voltage for different magnetic fields and tabulate the readings.
- 7. Measure the magnetic field using Gauss meter
- 8. From the data, calculate the Hall coefficient, carrier mobility and current density.

Table 1.1 Measurement of Hall coefficient

Current in the Hall effect setup = 2 mA

Current in the constant current power supply (A)	Magnetic field (H) (Gauss)	Hall Voltage (V _H) (volts)	Hall coefficient (R _H) cm ³ C ⁻¹

Observations and calculations

1. Thickness of the sample (t)	= cm
2. Resistivity of the sample (ρ)	$=$ $V C^{-1}$ s cm
3. Conductivity of the sample (σ)	$= \dots CV^{-1} s^{-1} cm^{-1}$

4. The hall coefficient of the sample $= R_H = \frac{V_H \cdot t}{IH} \times 10^8 = \dots \text{cm}^3 \text{ C}^{-1}$

5. The carrier density of the sample (n) $= \frac{1}{R_H q} = \dots \text{ carriers / cm}^{-3}$

6. The carrier mobility of the sample (μ) = $R_H \sigma$ = cm²V⁻¹ s⁻¹

Result

References

- 1. Ramsden and Edward Ramsden. Hall- effect Sensors. Amsterdam: Elsevier / Newnes. 2006
- 2. Fundamentals of solid state engineering by Manijeh Razeghi.
- 3. Materials handbook. Berlin:Springer, 2008.

2. Band Gap Determination using Post Office Box

Aim

To find the band gap of the material of the given thermistor using post office box.

Apparatus required

Thermistor, thermometer, post office box, power supply, galvanometer, insulating coil and glass beakers.

Principle

Wheatstone's Principle for balancing a network
$$\frac{P}{Q} = \frac{R}{S}$$

Of the four resistances, if three resistances are known and one is unknown, the unknown resistance can be calculated.

Formula

The band gap for semiconductors is given by
$$E_g = 2k \left(\frac{2.303 log_e R_T}{\frac{1}{T}} \right) eV$$

where
$$k = \text{Boltzmann constant} = 1.38 \times 10^{-23} \text{ J/K}$$

 $R_T = \text{Resistance at T K}$

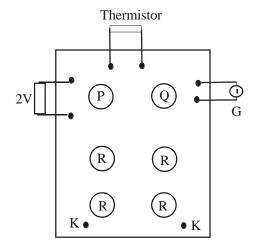


Fig.2.1Post Office Box - Circuit diagram

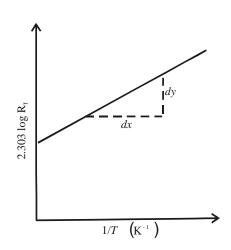


Fig.2.2 Model Graph

Procedure

- 1. The connections are given as in the Fig. 6.1(a).1. Ten ohm resistances are taken in P and Q.
- 2. Then the resistance in R is adjusted by pressing the tap key, until the deflection in the galvanometer crosses zero reading of the galvanometer, say from left to right.
- 3. After finding an approximate resistance for this, two resistances in R, which differ by 1 ohm, are to be found out such that the deflections in the galvanometer for these resistances will be on either side of zero reading of galvanometer.

- 4. We know $R_T = \frac{Q}{P} \times R = \frac{10}{10} \times R_1$ or $(R_1 + 1)$. This means that the resistance of the thermistor lies between R_1 and $(R_1 + 1)$. Then keeping the resistance in Q the same, the resistance in P is changed to 100 ohm.
- 5. Again two resistances, which differ by one ohm are found out such that the deflections in the galvanometer are on the either side of zero. Therefore the actual resistance of thermistor will be between $\frac{R_2}{10}$ and $\frac{R_2+1}{10}$.
- 6. Then the resistance in P is made 1000 ohms keeping same 10 ohms in Q. Again, two resistances R and (R+1) are found out such that the deflection in galvanometer changes its direction. Then the correct resistance = $R_T = \frac{10}{1000}(R)$ (or)R+1 = 0.01R (or) 0.01(R+1. Thus, the resistance of the thermistor is found out accurately to two decimals, at room temperature. The lower value may be assumed to be R_T (0.01R).
- 7. Then the thermistor is heated, by keeping it immersed in insulating oil. For every 10 K rise in temperature, the resistance of the thermistor is found out, (i.e) R_T 's are found out. The reading is entered in the tabular column.

Table 2.1To find the resistance of the thermistor at different temperatures

Temp. of thermistor T = t+273	$\frac{1}{T}$	Resistance in P	Resistance in Q	Resistance in R	Resistance of the thermistor $R_{T} = \frac{P}{Q}$ ×R	2.303 log ₁₀ R _T
K	K ⁻¹	ohm	ohm	Ohm	ohm	ohm

Observation

From graph, slope = (dy/dx) =

Graph

A graph is drawn between $\frac{1}{T}$ in X axis and 2.303 log R_T in Y axis where T is the temperature in K and R_T is the resistance of the thermistor at TK. The graph is shown in the Fig 2.2

Band gap (E_g)=2k × slope of the graph =
$$2k \times \left(\frac{dy}{dx}\right)$$

wherek = Boltzmann's constant.

Calculation

Band gap,
$$E_g = 2k(dy/dx) = \dots eV$$

Result

The approximate band gap value of given thermistor is $E_g = \dots eV$

- 1. B.L.Theraja, Basic Electronics Solid State, Fifth edition: S.Chand and company Ltd.
- 2. R.S.Sedha, A Textbook of Applied Electronics, Third edition: S.Chand and company Ltd.

3. To study V-I Characteristics of a Light Dependent Resistor (LDR)

Aim

To measure the photoconductive nature and the dark resistance of the given light dependent resistor (LDR) and to plot the characteristics of the LDR.

Apparatus required

LDR, Resistor (1 k Ω), ammeter (0 – 10 mA), voltmeter (0 – 10 V), light source, regulated power supply.

Principle

The photoconductive device is based on the decrease in the resistance of certain semiconductor materials when they are exposed to both infrared and visible radiation.

The photoconductivity is the result of carrier excitation due to light absorption and the figure of merit depends on the light absorption efficiency. The increase in conductivity is due to an increase in the number of mobile charge carriers in the material.

Formula

By ohm's law,
$$V = IR$$
 (or) $R = \frac{V}{I}$ ohm

where *R* is the resistance of the LDR (i.e.) the resistance when the LDR is closed. *V* and *I* represents the corresponding voltage and current respectively.

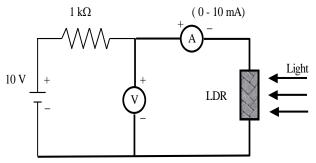


Fig. 3.1 Circuit diagram

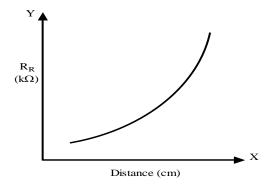


Fig.3.2 Model graph for photoconductivity of LDR

Procedure

- 1. The connections are given in as shown in Fig. 3.1
- 2. The light source is switched on and made to fall on the LDR.
- 3. The corresponding voltmeter and ammeter readings are noted.
- 4. The procedure is repeated by keeping the light source at different distances from the LDR.
- 5. A graph is plotted between resistance and distance of LDR from the light source.
- 6. The LDR is closed and the corresponding voltmeter and ammeter readings are noted. The value of the dark resistance can be calculated by Ohm's law.

Observation and calculation

Table 3.1 To determine the resistances of LDR at different distances

S.No	Distance	Voltmeter reading	Ammeter reading	R
5.110	(cm)	(V) volt	(I) mA	kΩ

Result

The characteristics of LDR were studied and plotted.

The dark resistance of the given LDR =ohm.

- 1. U.A Bakshi, A.P Godse, "Electronic Devices. Technical Publications", India, 2008.
- 2. Robert Diffenderfer, "Electronic Devices: Systems and Applications", Thomas Delmer Learning, USA, 2005.

4. Resistivity Determination for a Semiconductor Wafer using Four Probe Method

Aim

To determine the energy band gap of a semiconductor (Germanium) using four probe method.

Apparatus required

Probes arrangement (it should have four probes, coated with zinc at the tips). The probes should be equally spaced and must be in good electrical contact with the sample), Sample (Germanium or silicon crystal chip with non-conducting base), Oven (for the variation of temperature of the crystal from room temperature to about 200°C), A constant current generator (open circuit voltage about 20V, current range 0 to 10mA), Milli-voltmeter (range from 100mV to 3V), Power supply for oven Thermometer.

Formula

The energy band gapof semi-conductor is given by $E_g = 2k_B \frac{2.3026 \times log_{10} \rho}{\frac{1}{T}}$ in eV

where k_B is Boltzmann constant equal to 8.6×10^{-5} eV / kelvin pis the resistivity of the semi-conductor crystal given by

$$\rho = \frac{\rho_0}{f(W/S)} \text{ where } \rho_0 = \frac{V}{I} \times 2\pi s \text{ ; } \rho = \frac{V}{I} (0.213)$$

Here S is distance between probes

W is the thickness of semi-conducting crystal.

V and I are the voltage and current across and through the crystal chip.

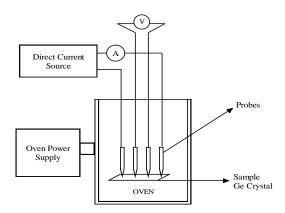


Fig 4.1 Four Probe experimental setup

Procedure

- 1. Connect one pair of probes to direct current source through milliammeter and other pair to millivoltmeter.
- 2. Switch on the constant current source and adjust current I, to a described value, say 2 mA
- 3. Connect the oven power supply and start heating.
- 4. Measure the inner probe voltage V, for various temperatures.

Table 4.1 To determine the resistivity of the semi-conductor for various temperatures

Current (I) = \dots mA

S.No.		erature	Voltage (V) Resistivity ρ (ohm. cm) 1/T (1)		1/T (10 ⁻³)	Log ₁₀ ρ
5.110.	in C	in K	(mV)	(K)		Log ₁₀ p

Observations

Graph

Plot a graph in $\left(\frac{10^3}{T}\right)$ and $\log_{10}\rho$ as shown in Fig.4.2. Find the slope of the curve $\frac{AB}{BC} = \frac{\log_{10}\rho}{10^3/T}$. So the energy band gap of semiconductor (Germanium) is given by $E_g = 2k \times \frac{2.3026 \times \log_{10}\rho}{1/T} = 2 \times 8.6 \times 10^{-5} \times 2.3026 \times \frac{AB}{BC} \times 1000 \text{ eV} = 0.396 \times \frac{AB}{BC} \text{ eV}$

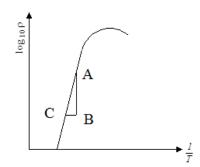


Fig. 4.2.Model Graph of $(10^3/T)$ and $log_{10}\rho$

Source of error and precautions

- 1. The resistivity of the material should be uniform in the area of measurement.
- 2. The surface on which the probes rest should be flat with no surface leakage.
- 3. The diameter of the contact between the metallic probes and the semiconductor crystal chip should be small compared to the distance between the probes.

Result

The energy band gap for given semiconductor is estimated to be $E_g = \dots eV$

- 1. Introduction To Semiconductor Materials And Devices, M.S Thyagi, Publisher John Wiley& Sons,2008
- 2. Physics of Semiconductor Devices, 3rd Ed, S.M.Sze, Kwok K .Ng, Publisher John Wiley& Sons,2008

5. Study of V-I and V-RCharacteristics of a Solar Cell

Aim

To study the V-I and V-R characteristics of a solar cell.

Apparatus Required

Solar cell, voltmeter, milliammeter, a dial type resistance box, keys, illuminating lamps, connecting wires etc.

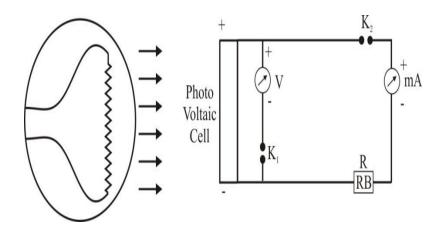


Fig.5.1 Schematic representation and circuit of Solar Cell

Procedure

A solar cell (photovoltaic cell) essentially consists of a p-n junction diode, in which electrons and holes are generated by the incident photons. When an external circuit is connected through the p-n junction device, a current passes through the circuit. Therefore, the device generates power when the electromagnetic radiation is incident on it. The schematic representation of a solar cell and the circuit connections are as shown in Fig. 5.1.

The voltmeter is connected in parallel with the given solar cell through a plug key. A milli-ammeter and a variable resistor are connected in series to the solar cell through a key as shown in the Fig.5.1.The solar cell can be irradiated by sun's radiation. Instead, it can also be irradiated by a filament bulb (60 W or 100 W). The resistance value is adjusted by a resistance box and the variation of V-I and V-R are plotted.

Table 5.1 V-I and V-R characteristics

Intensity	Resistance	Voltmeter Reading	Ammeter Reading
Maximum			

Table 5.2 V-I and V-R characteristics

Intensity	Resistance	Voltmeter Reading	Ammeter Reading
Minimum			

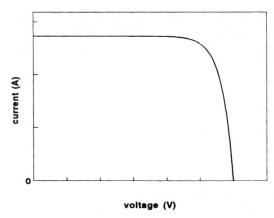


Fig. 5.2 Model Graph for V-I Characteristic

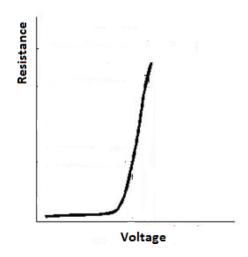


Fig. 5.3 Model Graph for V-R Characteristic

Result

The V-I and V-R Characteristics of the solar cell is studied.

- 1. Introduction To Semiconductor Materials And Devices, M.S Thyagi, Publisher John Wiley& Sons,2008
- 2. Physics of Semiconductor Devices, 3rd Ed, S.M.Sze, Kwok K .Ng, Publisher John Wiley& Sons,2008

6. Characteristic of PN junction diode under forward bias

Aim

To plot the characteristics curve of PN junction diode in Forward bias

Apparatus

A diode, DC voltage supplier, Bread board, 100Ω resistor, two multi-meters for measuring current and voltage and connecting wires

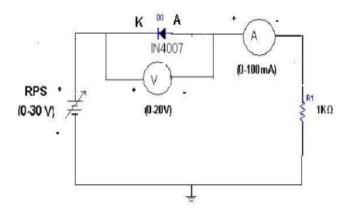


Fig. 6.1 PN Junction Forward Bias

Procedure:

For the forward bias of a P-N junction, P-type is connected to the positive terminal while the N-type is connected to the negative terminal of a battery. The potential at P-N junction can be varied with the help of potential divider. At some forward voltage (0.3 V for Ge and 0.7V for Si) the potential barrier is altogether eliminated and current starts flowing. This voltage is known as threshold voltage (V_{th}) or cut in voltage or knee voltage .It is practically same as barrier voltage V_{B} . For $V < V_{th}$, the current is negligible. As the forward applied voltage increase beyond threshold voltage, the forward current rises exponentially.

Table 1.1 V-I Characteristics of PN junction diode under forward bias

S. No	Forward Voltage (V _f) Volt	Forward Current (I _f)µA

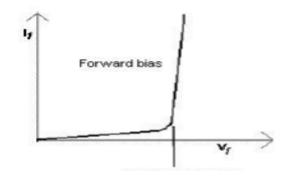


Fig. 6.2.Model graph for I-V Characteristics

Result

The PN junction characteristic is studied and curve is drawn.

- 1. B.L.Theraja, Basic Electronics Solid State, Fifth edition: S.Chand and company Ltd.
- 2. R.S.Sedha, A Textbook of Applied Electronics, Third edition: S.Chand and company Ltd.

7. To study V-I characteristics a photo cell

Aim

To study the illumination and V-I characteristic of photocell

Apparatus

Photo cell Apparatus and connecting wires.

Theory

Solid-state device with a photosensitive cathode that emits electrons when illuminated and an anode for collecting the emitted electrons. Illumination excites electrons, which are attracted to the anode, producing current proportional to the intensity of the illumination. In a photovoltaic cell, light is used to produce voltage. In a photoconductive cell, light is used to regulate the flow of current. Photocells are used in control systems, where interrupting a beam of light opens a circuit, actuating a relay that supplies power to a mechanism to bring about a desired operation, such as opening a door or setting off a burglar alarm.



Fig. 7.1 Experimental setup for Photo cell

Procedure

1. Illumination characteristics:

- 1. Connect the voltmeter and ammeter as per circuit diagram shown on the apparatus and set RL to minimum position.
- 2. A 100 W lamp is arranged over the solar panel such that light falls on it at normal. Initially, the lamp is placed at maximum distance and is switched on.
- 3. The voltage and the current is noted.
- 4. The intensity of the lamp is varied in steps by changing the distance of lamp. The readings are noted in Observation Table 7.1

2. Voltage – Current characteristics

- 1. The intensity of the lamp is kept at Medium (I_1) .
- 2. The voltage and current are switched on.
- 3. The load dial is set at 470 ohms. The voltage and current are noted down. The loads is set to 1 k, 2.2 k, 3.3 k, 4.7 k and 10 k ohms and at each setting the voltage and current are recorded. The readings are recorded in Observation Table 7.2

- 4. The intensity of the lamp is increased (I_2) . At this setting, step 3 is repeated.
- 5. The intensity of the lamp is set to its maximum (I_3) . At this setting step 3 is again repeated.

Table 7.1 Illumination characteristics


Intensity	At Minimum R_L ($k\Omega$)			
Variation (distance cm)	Voltage (V)	Current (μA)		

Table 7.2 V-I Characteristics of Photo cell

R _L		imum Intensity I ₁ ngest distance X ₁)		Intensity I ₂ distance X ₂)		n Intensity I3 st distance X3)
(kΩ)	Voltage (V)	Current (µA)	Voltage (V)	Current (µA)	Voltage (V)	Current (μA)

Graph

For V-I characteristics, a graph is plotted between the current and voltage

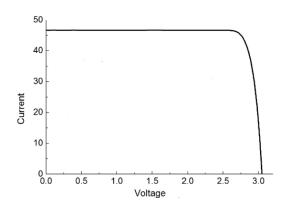

Result

Illumination and V-I characteristics of Photo cell has been studied.

- 1. BrijeshTripathi, Manoj Kumar, "Solar Energy: From Cells to Grid", CSFML Publications, India, 2018.
- 2. R.C. NevilleAngeleReinders, Pierre Verlinden, Wilfried van Sark, AlexandreFreundlich, "Photovoltaic Solar Energy: From Fundamentals to Applications", John Wiley & Sons Publication, UK, 2017.
- 3. R.C. Neville, "Solar Energy Conversion: The Solar Cell", Elsevier Science Publications, Netherlands, 1995.
- 4. Alan Fahrenbruch, Richard Bube, "Fundamentals of Solar Cells: Photovoltaic Solar Energy Conversion" Academic Press INC, UK, 1983.

8. Particle size determination using LASER

Aim

To determine the size of micro particles using laser.

Apparatus Required

Fine micro particles having nearly same size (say Lycopodium powder), a glass plate (say microscopic slide), diode laser, and a screen.

Principle

When laser is passed through a glass plate on which fine particles of nearly uniform size are spread, due to diffraction circular rings are observed. From the measurement of radii of the observed rings, we can calculate the size of the particles. Since for diffraction to occur size of the obstacle must be comparable with wavelength, only for extremely fine particles of micron or still lesser dimension, diffraction pattern can be obtained.

Diffraction is very often referred to as the bending of the waves around an obstacle. When a circular obstacle is illuminated by a coherent collimated beam such as laser light, due to diffraction circular rings are obtained as shown in the figure 8.1. If "r" is the radius of the first dark ring and "D" is the distance between the obstacle and screen on which the diffraction pattern is obtained, then.

$$\tan \theta = \frac{r}{D}$$

Since θ is very small in this experiment, $\tan \theta = \theta = \frac{r}{D}$

Formula

According to the theory, the diameter (2a) of the circular obstacle is given by

$$2a = \frac{1.22n\lambda D}{r_{-}}$$
 meters

where

 r_n =radius of the n^{th} order dark ring (m)

D = distance between the obstacle and the screen (m)

 λ =wavelength of the laser light (Å)

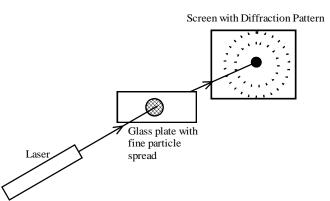


Fig. 8.1 Particle size determination using LASER

Procedure

- 1. Fine powder of particles is sprayed/spread on the glass plate.
- 2. Laser is held horizontally and the glass plate is inserted in its path.
- 3. Position of the glass plate is adjusted to get maximum contrast rings on the screen which is at a distance more than 0.5 m.
- 4. A white paper is placed on the screen and the positions of the dark rings are marked. The radii of different order dark rings (r_n) are measured using a scale.
- 5. The distance between the screen and the glass plate (D) is also measured. Using the given formula, the average diameter of the particles is calculated. $2a = \frac{1.22n\lambda D}{r}$
- 6. The experiment is repeated for different D values.

Table 8.1 Determination of particle size

S.No.	Distance (D)	Diffraction order (n)	Radius of dark ring (r _n)	Particle size (2a)
Unit	cm		cm	μm

Result

The average size of the particles measured using laser = $\dots \mu m$

- 1. Lasers: Fundamentals and Applications, Ajay Ghatak and K. Thyagarajan, Springer, 2010
- 2. Introduction to Optics, Frank L. Pedrotti , Leno M. Pedrotti, Leno S. Pedrotti, Cambridge University Press, 2017

9. Study of Attenuation and Propagation Characteristics of Optical Fiber Cable

9 (a) Attenuation in fibers

Aim

- 1. To determine the attenuation for the given optical fiber.
- 2. To measure the numerical aperture and hence the acceptance angle of the given fiber cables.

Apparatus Required

Fiber optic light source, optic power meter and fiber cables (1m and 5m), Numerical aperture measurement JIG, optical fiber cable with source, screen.

Principle

The propagation of light down dielectric waveguides bears some similarity to the propagation of microwaves down metal waveguides. If a beam of power P_i is launched into one end of an optical fiber and if P_f is the power remaining after a length L km has been traversed, then the attenuation is given by,

Attenuation =
$$\frac{10\log(\frac{P_i}{P_f})}{L} dB / km$$

Formula

Attenuation =
$$\frac{10\log\left(\frac{P_i}{P_f}\right)}{L} \text{ in dB / km}$$

where P_i = Power output for minimum light intensity (dB) P_f = Power output for maximum light intensity (dB)

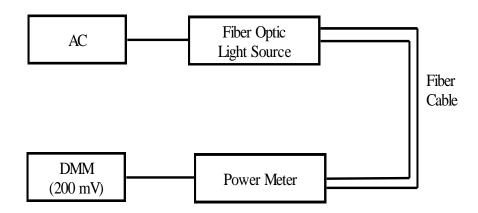


Fig. 9.1 Setup for loss measurement

Procedure

- 1. One end of the one metre fiber cable is connected to source and other end to the optical power metre.
- 2. Digital power meter is set to 200mV range (200 dB) and the power meter is switched on
- 3. The ac main of the optic source is switched on and the fiber patch cord knob in the source is set at one level (A).
- 4. The digital power meter reading is noted (P_i)
- 5. The procedure is repeated for 5m cable (P_f) .
- 6. The experiment is repeated for different source levels.

Table 9.1 Determination of Attenuation for optical fiber cables

$$L = 4 \text{ m} = 4 \times 10^{-3} \text{ km}$$

Source Level	Power output for 1m cable (P _i)	Power output for 5m cable (P _f)	Attenuation = $\frac{10\log(\frac{P_i}{P_f})}{L} dB / km$
Min			
Max			

Result

- 1. Attenuation at source level $A = \dots (dB/km)$
- 2. Attenuation at source level B = (dB/km)

9 (b) Numerical Aperture and Acceptance angle

Principle

Numerical aperture refers to the maximum angle at which the light incident on the fiber end is totally internally reflected and transmitted properly along the fiber. The cone formed by the rotation of this angle along the axis of the fiber is the cone of acceptance of the fiber.

Formula

Numerical aperture (NA)=
$$\frac{W}{\sqrt{4L^2+W^2}} = \sin\theta_{\text{max}}$$

Acceptance angle = $2 \theta_{max}$ (deg)

where L = distance of the screen from the fiber end in metre

W =diameter of the spot in metre.

 Θ = Half of acceptance angle (degree)

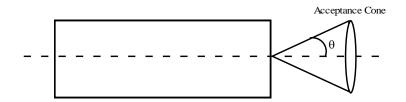


Fig. 9.2 Measurement of Numerical aperture

Procedure

- 1. One end of the 1 metre fiber cable is connected to the source and the other end to the NA jig.
- 2. The AC mains are plugged. Light must appear at the end of the fiber on the NA jig. The set knob in source is turned clockwise to set to a maximum output.
- 3. The white screen with the four concentric circles (10, 15, 20 and 25mm diameters) is held vertically at a suitable distance to make the red spot from the emitting fiber coincide with the 10mm circle.
- 4. The distance of the screen from the fiber end L is recorded and the diameter of the spot W is noted. The diameter of the circle can be accurately measured with the scale. The procedure is repeated for 15mm, 20mm and 25mm diameter circles.
- 5. The readings are tabulated and the mean numerical aperture and acceptance angle are determined.

Table 9.2 To find Numerical aperture and acceptance angle of optical fiber

Circle	Distance between source and screen (L) (mm)	Diameter of the spot W (mm)	$\mathbf{NA} = \frac{W}{\sqrt{4L^2 + W^2}}$	θ
5m				
1m				

Result

1. The Numerical aperture of fiber is measured as $1m = \dots \dots (No \ unit)$ $5m = \dots \dots (No \ unit)$ 2. The acceptance angle is calculated as $1m = \dots \dots (deg)$ $5m = \dots \dots (deg)$

- 1. Introduction to Fiber Optics, Ajay Ghatak and K. Thyagarajan, Cambridge University Press, 1998.
- 2. Fiber Optic Communications, Downing, James N., Clifton Park: Thomson Delmar Learning, 2004

10. Characteristic of PN junction diode under reverse bias

Aim

To plot the characteristics curve of PN junction diode in reverse bias

Apparatus

A diode, DC voltage supplier, Bread board, 100Ω resistor, 2 multimeter for measuring current and voltage and connecting wires.

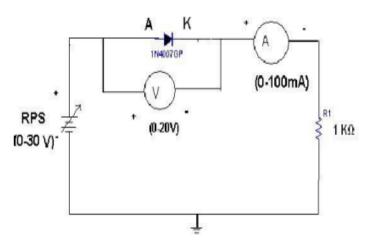


Fig. 10.1 Circuit diagram for PN Junction diode in Reverse Bias

Procedure

For the reverse bias of p-n junction, P-type is connected to the negative terminal while N-type is connected to the positive terminal of a battery. Under normal reverse voltage, a very little reverse current flows through a P-N junction. But when the reverse voltage is increased, a point is reached when the junction break down with sudden rise in reverse current. The critical value of the voltage is known as break down (VBR). The break down voltage is defined as the reverse voltage at which P-N junction breakdown with sudden rise in reverse current.

Table 10.1 Characteristics of PN junction diode under reverse bias

S. No	Reverse Voltage V _R Volt	Reverse Current I _R mA

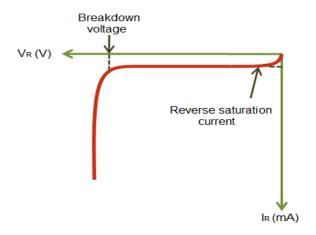


Fig. 10.2 Model graph for I-V Characteristics

Result

PN junction characteristic is studied and curve is drawn.

- 1. B.L.Theraja, Basic Electronics Solid State, Fifth edition: S.Chand and company Ltd.
- 2. R.S.Sedha, A Textbook of Applied Electronics, Third edition: S.Chand and company Ltd.

11. Determination of efficiency of Solar Cell

Aim

To explore solar cells as renewable energy sources and test their efficiency in converting solar radiation to electrical power.

Apparatus

Solar cell, Variable Resistor, Digital Multimeter (DMM), Electric motor, Desk lamp, Protractor Vernier Caliper.

Formula:

Efficiency of solar cell is measured as $\eta = [P_{max}/AI_o] \times 100$

where $P_{max} = Maximum power = I_{mp} x V_{mp}watt$

A = Area of the solar panel (7.2 cm x 4.5 cm)

 I_0 = Intensity of light = Power of the bulb/ $4\pi d^2$ W/m²

d = Distance between solar panel and bulb

Power of bulb = 75 watt

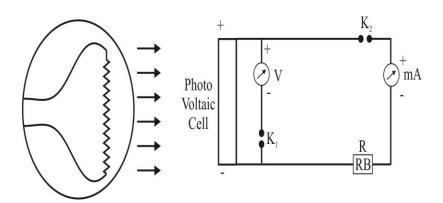


Fig. 11.1 Circuit diagram for efficiency measurement of solar cell

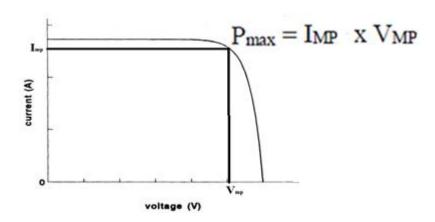


Fig. 11.2 Graph of V-I characteristics of solar cell

Procedure:

- 1. Measure the length and width of the solar cell with a vernier caliper and find its surface area
- 2. Place the desk lamp in front of solar panel
- 3. Measure the distance from solar cell to the desk lamp with the ruler. Adjust the distance to 0.05m and, turn on the desk lamp
- 4. Find the intensity of bulb as I_0 = Power of the bulb/ $4\pi d^2$
- 5. Connect the circuit as shown in the figure 11.1. A solar cell, variable resistance and a DMM (for measuring current) is connected in series. A second DMM (for measuring voltage) is connected to the solar cell in parallel.
- 6. Record the voltage and current values for different resistance values in table 11.1
- 7. Adjust the distance from solar cell to the desk lamp as 0.10cm and record the voltage and current for different resistance
- 8. Plot the graph with voltage along X-axis and current along Y-axis and find maximum power as $P_{max} = I_{mp} \ x \ V_{mp}$ watt

Table 11.1 V-I characteristicswhen d = 5 cm

Intensity	Resistance	VoltmeterReading	AmmeterReading
Maximum			

Table 11.2 V-I characteristics when distance d = 10 cm

Resistance	VoltmeterReading	AmmeterReading
	Resistance	Resistance VoltmeterReading

Observation

Intensity of light I_0 when d = 5 cm $= \dots W / m^2$ Intensity of light Io when d = 10 cm $= \dots W/m^2$ When d = 5cm, Power maximum $= \dots$ When d = 10cm, Power maximum $= \dots$

Result:

Efficiency of solar cell is calculated as, $\eta \; (\text{minimum intensity}) = \dots \dots \%$ $\eta \; (\text{maximum intensity}) = \dots \dots \%$

- 1. Introduction To Semiconductor Materials And Devices, M.S Thyagi, Publisher John Wiley& Sons,2008
- 2. Physics of Semiconductor Devices, 3rd Ed, S.M.Sze, Kwok K .Ng, Publisher John Wiley& Sons,2008

12. Calculation of Lattice Cell Parameters – X-ray Diffraction

Aim

Tocalculate the lattice cell parameters from the powder X-ray diffraction data.

Apparatus required

Powder X-ray diffraction diagram

Principle

Bragg's law is the theoretical basis for X-ray diffraction.

$$(\sin^2 \theta)_{hkl} = (\lambda^2 / 4a^2) (h^2 + k^2 + l^2)$$

Each of the Miller indices can take values 0, 1, 2, 3...... Thus, the factor $(h^2 + k^2 + l^2)$ takes the values given in Table 12.1

Formula

For a cubic crystal

$$\frac{1}{d^2} = \frac{(h^2 + k^2 + l^2)}{a^2}$$

For a tetragonal crystal

$$\frac{1}{d^2} = \left\{ \frac{(h^2 + k^2)}{a^2} + \frac{l^2}{c^2} \right\}$$

For a orthorhombic crystal

$$\frac{1}{d^2} = \left(\frac{h^2}{a^2}\right) + \left(\frac{k^2}{b^2}\right) + \left(\frac{l^2}{c^2}\right)$$

The lattice parameter and interplanar distance are given for a cubic crystal as,

$$a = \frac{\lambda}{2\sin\theta} \sqrt{h^2 + k^2 + l^2} \mathring{A}$$

$$d = \frac{a}{\sqrt{h^2 + k^2 + l^2}} \mathring{A}$$

where,

a = Lattice parameter

d = Interplaner distance

 λ = Wavelength of the CuK α radiation (1.5405)

h, k, l = Miller integers

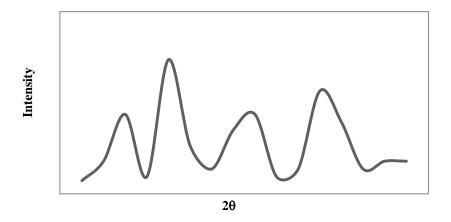


Fig. 12.1. Model X-ray diffraction pattern

The problem of indexing lies in fixing the correct value of a by inspection of the $\sin^2\theta$ values.

Procedure:

From the 20 values on a powder photograph, the θ values are obtained. The $\sin^2\theta$ values are tabulated. From that the values of $1 \times \frac{\sin^2\theta}{\sin^2\theta_{min}}$, $2 \times \frac{\sin^2\theta}{\sin^2\theta_{min}}$, $3 \times \frac{\sin^2\theta}{\sin^2\theta_{min}}$ are determined and are tabulated. The values of $3 \times \frac{\sin^2\theta}{\sin^2\theta_{min}}$ are rounded to the nearest integer. This gives the value of $h^2 + k^2 + l$. From these the values of h, k, l are determined from the Table 12.1

From the h,k,l values, the lattice parameters are calculated using the relation

$$a = \frac{\lambda}{2\sin\theta} \sqrt{h^2 + k^2 + l^2} \quad \mathring{A}$$
$$d = \frac{a}{\sqrt{h^2 + k^2 + l^2}} \quad \mathring{A}$$

Table 12.1 Value of $h^2 + k^2 + l^2$ for different planes

h, k, l	$h^2 + k^2 + l^2$	h, k, l	$h^2 + k^2 + l^2$
100	1	300	9
110	2	310	10
111	3	311	11
200	4	322	12
210	5	320	13
211	6	321	14
220	8	400	16
221	9	410	17

Table 12.2 Lattice parameter (a) and Interplanar distance (d) for different planes

S. No	2θ	$sin^2\theta$	$1 \times \frac{\sin^2 \theta}{\sin^2 \theta_{\min}}$	$2 \times \frac{\sin^2 \theta}{\sin^2 \theta_{\min}}$	$3 \times \frac{\sin^2 \theta}{\sin^2 \theta_{\min}}$	$h^2 + k^2 + l^2$	hkl	a Å	d Å

Table 3 Lattice determination

Lattice type	Rule for reflection to be observed
Primitive P	None
Body centered I	hkl: h + k + l = 2 n
Face centered F	hkl: h, k, l either all odd or all even

Depending on the nature of the h,k,l values the lattice type can be determined.

Result:

The lattice parameters are calculated theoretically from the powder x-ray diffraction pattern.

- 1. Dann, S.E. Reactions and Characterization of SOLIDS. Royal Society of Chemistry, USA (2002).
- 2. Skoog, D.A.; Holler, F.J.; Crouch, S.R. Principles of Instrumental Analysis. SixthEdition, Thomson Brooks/Cole, USA (2007).

13. Mini Project – Concept based Demonstration

Aim:

To construct the working model based on principles of physics, the opportunity to develop a range of skills and knowledge already learnt to an unseen problem.

Objectives:

On successful completion of the mini project, the student will have developed skills in the following areas:

- 1. Design of experiments.
- 2. Experimental or computational techniques.
- 3. Searching the physical and related literature.
- 4. Communication of results in an oral presentation and in a report.
- 5. Working as part of a team.
- 6. Assessment of team members.

Assessment and Evaluation:

- 1. Each Class should have at least **eight groups**. Each group should have a **minimum of 5 members** or above and Maximum of **9 members**.
- 2. Mini Project should be a **working model**. One page write-up about the project should be submitted as per the template provided by the class subject teacher.
- 3. **Department of Physics & Nanotechnology** will be organizing an event **TechKnow** to showcase these Mini Projects. All groups should present the working model along with the poster at the TechKnow.
- 4. Expert Committee will evaluate and select the best project from each class.
- 5. Certificates will be awarded for the best project during the event TechKnow.
- 6. Marks for the project will be awarded under the following criteria.

S. No	Criteria	Marks
1.	Working model / Design	5
2.	Idea/ Concept / Novelty	5
3.	Presentation / Viva	5
4.	Usefulness / Application	5
	Total	20